

OXYPLUS™ INVESTMENT & LIQUID



Ideal for Casting Low Expanding Patterns

Oxyplus investment is a versatile oxyphosphate dental investment and liquid combination. Formulated to provide easy removal and smooth castings. Oxyplus investment is ideal for making refractory models for waxing and casting precious and nonprecious, high temperature alloys for partial denture frames.

Oxyplus investment is designed for use with both Oxyplus liquid and distilled water, allowing for flexibility in application. With its longer setting time, it ensures accuracy and reliability in dental casting applications.

Color: Off white

Ratios

When making the refractory model, Oxyplus liquid is recommended at a diluted ratio of 3 parts deionized or distilled water to 1 part Oxyplus liquid by volume. Oxyplus liquid ratios will need to be slightly adjusted based on your geographical location and/or climate conditions.

Please note:

- More Oxyplus liquid = more expansion for looser fitted castings
- Less Oxyplus liquid = less expansion for tighter fitted castings

When investing the outer mold, Oxyplus investment is mixed with deionized or distilled water only. Do not use Oxyplus liquid.

Duplicating the Gypsum Master Model

1. Soak the master model in water for 30 minutes (or 2 minutes under vacuum) to remove any air from the gypsum model.
2. Blow off any excess water from the model and place the model into a duplicating flask.
3. Pour Uniflex™, Multiflex™, Perfex®, or Instaflex™ Concentrate duplicating material from your duplicating unit.
4. Once filled, place the duplicating flask into a chilled water bath, covering only the bottom of the flask and not the sides.
5. As the duplicating material starts to cool and set up, it will pull toward the cold water, which will cause it to pull tightly over the master model for an exact duplication.
6. Let set for 30 minutes, then remove the master model and allow the duplicating material to bench set to reach room temperature before pouring in Oxyplus investment.

Making the Oxyplus Refractory Model

	Mechanical Mix	Hand Mix
Powder	200g	200g
Diluted Oxyplus Liquid (3 to 1 Ratio)	25 ml (18.75 ml distilled water to 6.25 ml Oxyplus liquid)	25 ml (18.75 ml distilled water to 6.25 ml Oxyplus liquid)
Mix Time	60 seconds under vacuum	90 seconds

1. Place impression on vibrator and fill the impression with the mixture while vibrating.
2. Remove from vibrator and allow to harden for 30 minutes.
3. Remove the model, trim and dry for 45 minutes at 375-425 °F (190-218 °C).
4. Follow the application instructions for an R&R Rosin Dip and dip the model.



RANSOM & RANDOLPH

3535 Briarfield Boulevard | Maumee, OH 43537 USA
800.800.7496 | 419.865.9497 | 419.865.9997 (FAX)
www.ransom-randolph.com



Made in the USA

Ransom & Randolph GmbH
Leipziger Straße 40 | 04571 Rötha Germany
+49 342 06373999

Investing with Innovation™

OXYPLUS™ INVESTMENT & LIQUID


Investing the Model (Outer Mold)

	Mechanical Mix	Hand Mix
Powder	600g	600g
Deionized or Distilled Water	70 ml	73 ml
Mix Time	60 seconds under vacuum	90 seconds

1. Place the ring on the vibrator and pour the mixed investment into the ring until the wax-up is covered.
2. Vibrate the mix back out of the ring until there is a bubble-free thin coating covering the wax-up.
3. Once a bubble-free thin coating is achieved, vibrate the mix back into the ring until filled.
4. Remove immediately from vibrator.
5. Bench set for 60 minutes.

Burnout

From	To	In	Hold	Raise To	In	Hold
70 °F (21 °C)	550 °F (290 °C)	20 min	30 min	1800 °F (982 °C)	60 min	60 min

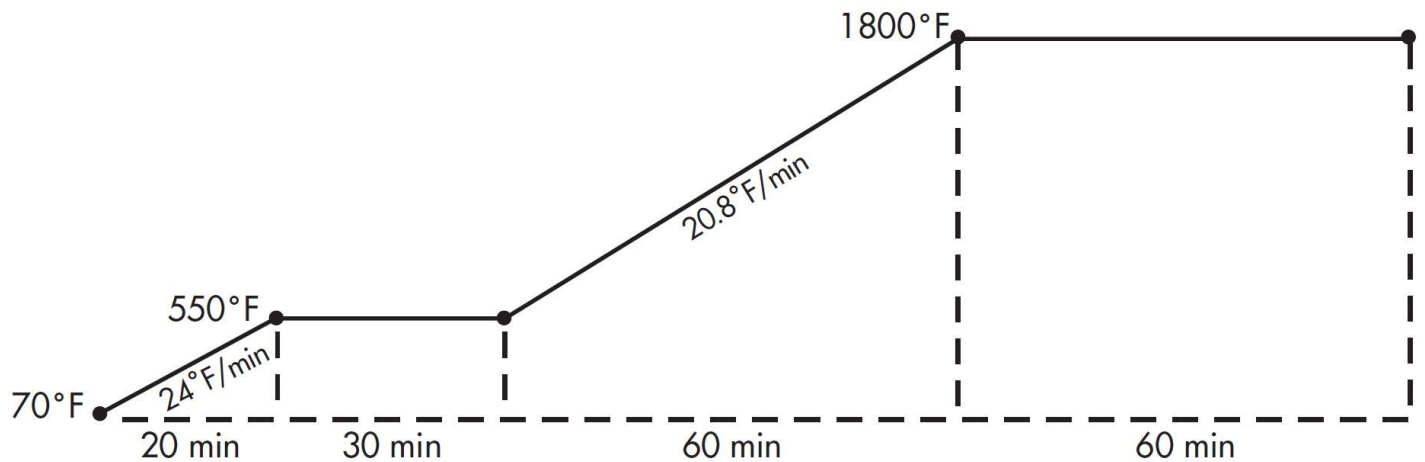

Technical Tips

For additional information and recommendations, refer to the Dental Investment Defect Analysis Guide Technical Tip available for download at ransom-randolph.com/technical-tips.

Storage

Keep container closed when not in use. Store in a cool dry area.

Safety

North America: Danger. Contains crystalline silica. May cause cancer by inhalation. Causes damage to lungs through prolonged or repeated exposure by inhalation. See SDS for more information.



RANSOM & RANDOLPH

3535 Briarfield Boulevard | Maumee, OH 43537 USA
 800.800.7496 | 419.865.9497 | 419.865.9997 (FAX)
www.ransom-randolph.com



Ransom & Randolph GmbH
 Leipziger Straße 40 | 04571 Rötha Germany
 +49 342 06373999

Investing with Innovation™

OXYPLUS™ INVESTMENT & LIQUID

EU: Danger. Contains respirable crystalline silica. Causes damage to lungs through prolonged or repeated exposure. See SDS for more information.

Availability

Oxyplus investment

Package Size	Item Number
400 gram pouch (quantity of 30)	400040
40 LB (18.2 KG) pail	400043

Oxyplus Liquid

Package Size	Item Number
1 GAL (3.8 L) jug	400047

To order, visit shop.ransomrandolph.com or contact an authorized distributor. For a list of authorized distributors near you, visit ransom-randolph.com/distributors or contact your region's R&R Customer Service Team.

As the conditions or methods of use are beyond our control, we do not assume any responsibility and expressly disclaim any liability for any use of this product. Information contained herein is believed to be true and accurate but all statements or suggestions are made without warranty, expressed or implied, regarding accuracy of the information, the hazards connected with the use of the material or the results to be obtained from the use thereof. Compliance with all applicable federal, state, and local regulations remains the responsibility of the user. All potential liability related to the sale and use of this product is limited to the cost of the particular goods sold in their respective transactions.



RANSOM & RANDOLPH

3535 Briarfield Boulevard | Maumee, OH 43537 USA
800.800.7496 | 419.865.9497 | 419.865.9997 (FAX)
www.ransom-randolph.com



Made in the USA

Issue Date: July 26, 2024

Ransom & Randolph GmbH
Leipziger Straße 40 | 04571 Rötha Germany
+49 342 06373999

Investing with Innovation™